A 506

de	Gebrauchsanweisung Wagen
da	Brugsanvisning Vogn
en	Operating instructions: Mobile unit
es	Instrucciones de manejo Carros
fi	Käyttöohje Vaunu
fr	Mode d'emploi Chariot
it	Istruzioni d'uso Carrelli
nl	Gebruiksaanwijzing Wagen
no	Bruksanvisning vogn
pt	Instruções de utilização Carro
sv	Bruksanvisning Vagn

de	4
da	15
en	26
es	37
fi	48
fr	59
it	70
nl	81
no	92
pt	103
sv	114

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Notes about these instructions	26
Intended use	27
Queries and technical problems.....	27
Items supplied	28
Load carriers	28
Disposal of the packing material	28
Optional accessories.....	29
Warning and Safety instructions	30
Medical devices.....	30
Laboratory glassware and utensils.....	30
Areas of application	31
To check before loading the machine and before starting a programme	31
Programme selection	31
Medical devices.....	32
Loading mesh trays	32
Reprocessing hollow items	33
Laboratory glassware and utensils.....	34
Watch glasses	34
Small items	34
Microtiter plates and blood analysis plates.....	35

Important warnings

 Information which is important for safety is highlighted in a thick framed box with a warning symbol. This alerts you to the potential danger of injury to people or damage to property.

Read these warning notes carefully and observe the procedural instructions and codes of practice they describe.

Notes

Information of particular importance that must be observed is highlighted in a thick framed box.

Additional information and comments

Additional information and comments are contained in a simple frame.

Operating steps

Operating steps are indicated by a black square bullet point.

Example:

- Select an option using the arrow buttons and save your choice with *OK*.

Display

Certain functions are shown in display messages using the same font as that used for the function itself in the display.

Example:

Settings  menu.

This mobile unit can be used to reprocess machine-reprocessible medical devices, laboratory glassware and utensils in a Miele washer-disinfector. Follow the operating instructions for your washer-disinfector as well as the instructions of the medical device or glassware and utensil manufacturer on how to reprocess their items by machine.

The A 506 mobile unit can accommodate up to 12 DIN mesh trays or can be loaded with a variety of mesh trays and inserts for reprocessing medical devices, laboratory glassware and utensils.

The mobile unit can be used in the following washer-disinfectors:

- PLW 8615
- PLW 8616
- PLW 8617
- PWD 8626
- PWD 8628

The washer-disinfector is generally referred to as “the machine” in these operating instructions. Reprocessible medical devices, laboratory glassware and utensils are referred to as “items” if they are not more closely defined.

Queries and technical problems

If you have any queries or technical problems that you would like to discuss please contact Miele. Contact details can be found at the back of the Operating Instructions for your washer-disinfector or go to www.miele-professional.com.

Load carriers

A 506 mobile unit, height 557 mm, width 650 mm, depth 790 mm

A 506 *	Loading dimensions [height x width x depth in mm]
Level 1	112.5 x 585 x 780
Level 2	112.5 x 595 x 780
Level 3	112.5 x 595 x 780
Level 4	114.0 x 595 x 780

* Levels are numbered from the bottom up

Disposal of the packing material

The packaging is designed to protect against transportation damage. The packaging materials used are selected from materials which are environmentally friendly for disposal and should be recycled.

Recycling the packaging reduces the use of raw materials in the manufacturing process and also reduces the amount of waste in landfill sites.

Optional accessories

The following and other accessories are available to order from Miele:

- A 810 Cover net, metal frame with replaceable silicone net A 811, 428 x 220 mm
- E 142, DIN mesh tray with a 5 mm grid and 2 hinged handles, height 45/55 mm, width, 255 mm, depth 480 mm
- E 146, mesh tray with lid and hinged handles, height 55 mm, width 150 mm, depth 255 mm

Medical devices

- E 473/2 mesh tray with lid for small items, height 85 mm, width 60 mm, depth 60 mm
- E 476, 50 holders for Ø 4 to 8 mm micro instruments, for use in mesh trays with a 5 mm grid
- E 477, 20 holders for micro-instruments, for use in mesh trays with a 5 mm grid
- E 479, 50 holders for Ø 4 mm micro instruments, for use in mesh trays with a 5 mm grid

Laboratory glass-ware and utensils

- E 402, insert for 44 watch glasses Ø 80 - 125 mm, spacing approx. 15 mm, height 53 mm, width 200 mm, depth 445 mm.
- E 403, insert for 105 watch glasses Ø 50-60 mm, spacing approx. 9 mm, height 35 mm, width 200 mm, depth 445 mm.
- E 494, insert for 5 loose microtiter plates, height 35 mm, width 205 mm, depth 440 mm.

en - Warning and Safety instructions

To avoid the risk of accidents and damage to this mobile unit please read these instructions carefully before using it for the first time.
Keep these instructions in a safe place where they are accessible to users at all times.

 Please also read the operating instructions for your washer-disinfector and pay particular attention to the Warning and Safety instructions.

- ▶ The mobile unit is approved solely for the applications specified in the “Areas of application” chapter of these operating instructions. Applications that go beyond those listed and use of the mobile unit with additional components from the Miele product range are described in the operating instructions for the components or are to be agreed with Miele.
- ▶ New mobile units and baskets must be washed in the washer-disinfector without a wash load prior to first use.
- ▶ Inspect all mobile units, baskets, modules and inserts daily as described in the “Maintenance” section in the Operating instructions for the washer-disinfector.
- ▶ Broken glass can result in serious injury when loading or unloading. Damaged glass items must not be reprocessed in the machine.

Medical devices

- ▶ The reprocessing programme used should include thermal disinfection, and the chemical agents used must be suitable for the material quality of the wash load.
- ▶ The last rinse, as a minimum, must be carried out with DI water.
- ▶ Visual checks should be undertaken to check the cleaning result. Additional protein analysis checks should also be carried out e.g. weekly, on a random sampling basis.

Laboratory glassware and utensils

- ▶ The reprocessing result must be checked as necessary using a suitable test method and not just by a visual check.

Miele cannot be held liable for damage caused by non-compliance with these Warning and Safety Instructions.

For validation purposes please make sure that you follow the loading instructions given on the template.

To check before loading the machine and before starting a programme

- Are the spray arms able to rotate freely?
- Check spray arm jets for blockages.
- Check that the mobile unit is correctly docked to the water path in the chamber.

Programme selection

When selecting a programme manually on the machine, the relevant mobile unit must be selected first. Pictograms representing the mobile units are shown on the display.

The pictogram can be found in the bottom right on the front of the mobile unit, next to the RFID tag.


For the A 506 mobile unit, the pictogram below is displayed.

Once you have selected the pictogram, the list of assigned reprocessing programmes is displayed.


Mobile unit	Pictogram
A 506	

Medical devices

Loading mesh trays

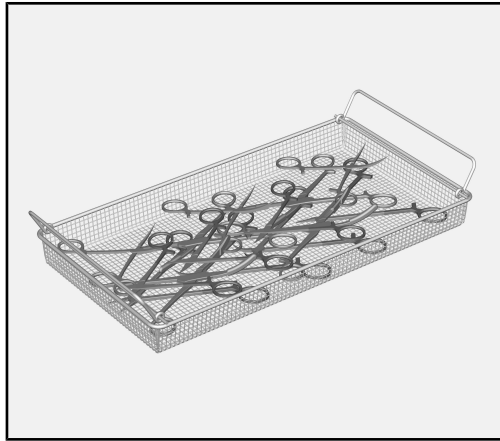
- Hinged instruments must be opened to minimise the risk of surfaces not being accessed during reprocessing.

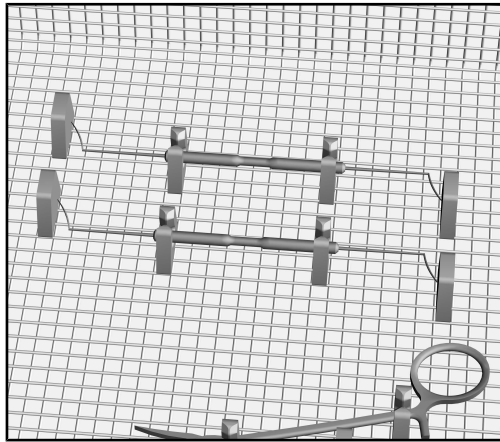- To prevent the risk of light instruments and other small items being moved about during reprocessing, place them in closable inserts such as the E 473/2 or secure them in silicone holders such as the E 476, E 477 and E 479. The E 473/2 insert is supplied with its own operating instructions.

Reprocessing hollow items

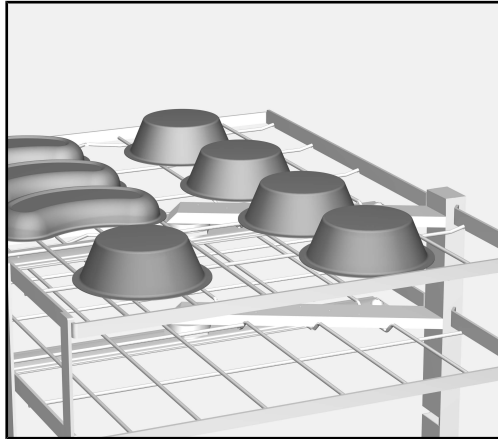- Hollow items such as beakers and kidney dishes should be loaded with the open side facing downwards. Lightweight items should be secured with a cover net, e.g. the A 810.
- Special inserts should be used if necessary for reprocessing hollow items. These inserts can be purchased from Miele.

The hollow items must not hinder the reprocessing of other items and must be secured so that they cannot be flipped over by the force of water during reprocessing.

If hollow items are placed with the open side facing upwards, suds solution can collect in them. This can spill and impair reprocessing results.

Laboratory glassware and utensils

Items with wide necks, e.g. glass beakers, wide-necked Erlenmeyer flasks and petri dishes, or cylindrical items, e.g. test tubes, can be cleaned and rinsed inside and out using rotating spray arms.

- Load wide-necked glassware with the open side facing downwards.

The laboratory glassware must not hinder the reprocessing of other items and must be secured so that it cannot be flipped over by the force of water during reprocessing. Secure light items with a cover net such as the A 810 or reprocess them in mesh trays with lids.

If laboratory glassware is placed with the open side facing upwards, suds solution can collect in it. This can spill and impair reprocessing results.

- Special inserts should be used if necessary for reprocessing laboratory glassware and utensils. Various inserts for use with different items can be purchased from Miele. Some examples are given below.

Watch glasses

Special inserts are available for reprocessing watch glasses, e.g. the E 403 for up to 105 watch glasses.

Small items

Reprocess small items such as stoppers, lids, spatulas, magnetic stirring rods, etc. in mesh trays with lids, e.g. the E 146. Alternatively a mesh tray with a sufficiently large cover net can be used.

Microtiter plates and blood analysis plates

The E 494 insert is designed for the machine reprocessing of microtiter plates and blood analysis plates.

Plastic microtiter plates and blood analysis plates:

For microtiter plates / blood analysis plates made of polycarbonate, polyacrylic, polystyrene, etc., please note the following:


- Use cleaning agents with hypochlorite or hydrogen peroxide oxidising properties.

⚠ Do not use cleaning agents and defoamers containing surfactants and/or silicone oils.

- Use neutralising agents with a citric acid base.
- Select a cleaning programme with 1 - 2 cold pre-washes.
- The temperature must not exceed 55 °C in the washing and rinsing blocks.
- After the final rinse, manually tip any remaining water out of the cavities or wells before drying.
- The temperature when drying in the washer-disinfector must not exceed 65 °C.
- The temperature when drying in the drying cabinet must not exceed 55 °C.

The instructions of the manufacturer of the microtiter plates or the blood analysis plates should also be observed.

Arrange microtiter plates in the E 494 insert

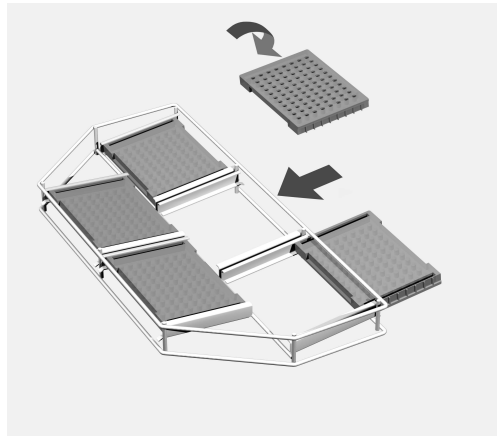- Push the microtiter plates into the insert holders with the open sides facing downwards.
- Up to two inserts can be used per level. Do not stack the inserts on top of each other.

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