

ProCare

Dent 40

Rinsing agent

Use:

Liquid concentrate. For the automated reprocessing of dental instruments including hand pieces in Miele Professional washer disinfectors.

Quality features:

ProCare Dent 40 complies with European guidelines for medical devices. **CE MD**

If a serious incident occurs with the product, report it to the manufacturer and the relevant national authority.

Product properties:


- Improves and accelerates the drying of dental instruments in machine-based reprocessing methods.
- Neutralises alkaline carryover and alkalinity.

Use and dispensing:

Miele *ProCare Dent 40* is formulated for use in processes in Miele Professional washer disinfectors.

- Suitable dosing dispensers should be used.
- Application concentration: 0.3 ml/l (0.03 %) via door dispenser for rinsing agents
- Program selection: Vario TD

Program sequence:

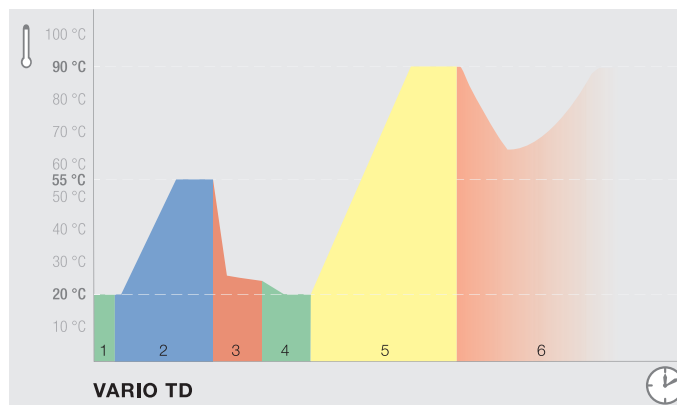

Vario TD process:

- 1 Pre-wash with cold water to remove coarse soiling (e.g. blood)
- 2 Main wash
- 3 Neutralisation
- 4 Intermediate rinse with cold water
- 5 Thermal disinfection at 93 °C with 5-minute holding time in the final rinse with *ProCare Dent 40* (dispensing temperature 70 °C)
- 6 Drying

Miele recommends the use of fully demineralised water in the final rinse in order to protect the load and the washer-disinfector against stains and deposits and to protect anodised aluminium.

Please note!

- Do not mix with other products.
- Call Miele Service when switching products.
- Reprocessing should comply with all ordinances pursuant to medical device regulations and should be performed with appropriate validated processes.
- Observe recommendations on reprocessing issued by medical device manufacturers with respect to the provisions of CAN/CSA-Z17664 (ISO 17664).
- Operating instructions on Miele Professional washer disinfectors must be observed.
- For commercial applications only.

ProCare

Dent 40

Rinsing agent


Technical data:

- Density (at 20 °C): approx. 1.1 g/cm³
- pH-value (determined using fully demineralised water at 20 °C) 0.3 ml/l: approx. 3.8

Ingredients according to Regulation (EC) No 648/2004 on detergents:

- < 5 % phosphonates
- 15 % to 30 % non-ionic tensides
- Also preservatives (methylchloroisothiazolinone / methylisothiazolinone)

Storage information:

- Storage temperature: 0 °C to 30 °C
- Can be used for 3 years under proper storage conditions
- Use-by date: See label next to hourglass symbol 

Warning and Safety notes:

For safety and environmental information, cf. safety data sheets. These are available from the Miele websites listed below.

Dispose of containers only when empty and sealed. Disposal of product residue: refer to safety data sheet.

All information contained herein is current and correct to our knowledge. This does not however release users from obligations to perform their own trials and experiments. This cannot be deemed as a legal warranty of product characteristics or properties.

Packaging unit:

Designation	Packaging unit	Mat. no.
ProCare Dent 40	1 l bottle	12 187 410

Contact information:**Germany**

Miele & Cie. KG
Carl-Miele-Straße 29,
33332 Gütersloh
Tel: 49 (0) 8002244644
vertrieb@miele-professional.de
www.miele.de/professional

Canada

Miele Limited
161 Four Valley Drive
Vaughan, ON L4K 4V8
Tel: 1-888-325-3957
email: professional@miele.ca
www.mieleprofessional.ca

**Germany**

Chemische Fabrik
Dr. Weigert GmbH & Co. KG
Mühlenhagen 85
20539 Hamburg
www.drweigert.de