A 207

cs	Návod k obsluze Injektorový vozík
el	Οδηγίες χρήσης ειδικό συρόμενο βαγονέτο με εγχυτήρες
en	Operating instructions Injector unit
hu	Használati utasítás Injektoros kocsi
ja	取扱説明書 インジェクターユニット
ko	사용 설명서 인젝터 장치
pl	Instrukcja użytkowania Wózek iniekcyjny
ro	Instrucțiuni de utilizare Cărucior cu injector
ru	Инструкция по эксплуатации Инжекторная тележка
sk	Návod na použitie Injektorový vozík
tr	Kullanım kılavuzu Mobil enjektör ünitesi
zh	操作说明书 注射装置

cs	4
el	16
en	28
hu	40
ja	52
ko	64
pl	76
ro	88
ru	100
sk	112
tr	124
zh	136

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Notes about these instructions	28
Queries and technical problems	28
Intended use	29
Items supplied	30
Scope of delivery A 207	30
Disposal of the packing material	31
Optional components	31
Warnings and safety notes	32
Areas of application	33
Programme selection	33
To check before loading the machine and before starting a programme	33
Loading ophthalmic instruments	34
Loading mesh trays without injector manifolds	35
A 314 mesh tray	37
Testing lumened instruments within the framework of quality assurance tests	37
Test point for measuring water pressure	38

Important warnings

 Information which is important for safety is highlighted in a thick framed box with a warning symbol. This alerts you to the potential danger of injury to people or damage to property. Read these warning notes carefully and observe the procedural instructions and codes of practice they describe.

Notes

Information of particular importance that must be observed is highlighted in a thick framed box.

Additional information and comments

Additional information and comments are contained in a simple frame.

Operating steps

Operating steps are indicated by a black square bullet point.

Example:

- Select an option using the arrow buttons and save your choice with *OK*.

Display

Certain functions are shown in display messages using the same font as that used for the function itself in the display.

Example:

Settings  menu.

Queries and technical problems

In the event of queries or technical problems, please contact Miele. Contact details can be found at the end of the operating instructions for your cleaning machine or at www.miele.com/professional.

This mobile unit can be used to reprocess machine-reprocessible medical devices in a Miele washer-disinfector. Follow the operating instructions for your washer-disinfector as well as the instructions issued by the medical device manufacturer on how to reprocess their items using a machine.

The A 207 mobile injector unit is intended for reprocessing ophthalmic instruments only and may only be used in the OphthaTrays A207 programme.

This mobile injector unit has 3 levels and 2 spray arms.

The top level is fitted with an injector manifold with silicone hoses and Luer Lock connectors. It can be used for connecting mesh trays, baskets, cassettes, bowls, etc with integrated injector manifolds for ophthalmic sets.

The two lower levels are designed to take inserts and mesh trays for reprocessing instruments without lumen.

Narrow lumened instruments (Sauter cannulae or rinsing and suction hand pieces, for example) require the suds solution to be filtered. The injector manifolds must be fitted with a Miele A 800 tubular filter for this purpose. The tubular filter is supplied with its own operating instructions.

Special reprocessing programmes are required for reprocessing instruments from the above areas of application by machine. These programmes are installed ex-works or can be installed, created or activated by the Miele Customer Service Department on the following machines:

- | | |
|---------------|--------------|
| - PWD 8682 | - PG 8582 |
| - PWD 8682 CD | - PG 8582 CD |
| - PWD 8692 | - PG 8592 |
| - DS 5010 | |
| - DS 5011 | |
| - DS 5020 | |

Disposable items must not be reprocessed.

The washer-disinfector is generally referred to as “the machine” in these operating instructions. Reprocessible medical devices are referred to as “load items” if they are not more closely defined.

Mesh trays, baskets, cassettes, bowls etc and other preparation containers are referred to throughout these instructions as mesh trays.

Scope of delivery A 207

- Mobile injector unit A 207
 - Height 438 mm, width 530 mm, depth 542 mm
 - 3 levels
 - 1 nozzle manifold with 11 connections for irrigation connectors
 - Spray arms with magnets for spray arm monitoring
- A 800, reusable tubular filter with a filtration grade of ≥ 0.07 mm, 405 mm long, $\varnothing$ 22.5 mm
- 2 x A 311, insert with 4 compartments for mesh bowls or trays, spacing 50 mm, Height 142 mm, Width 197 mm, Depth 436 mm
- 3 x E 362 blind stoppers
- 8 silicone hoses with Luer Lock ♂ adapters for Luer Lock ♀, length 160 mm

Disposal of the packing material

The packaging is designed to protect against transportation damage. The packaging materials used are selected from materials which are environmentally friendly for disposal and should be recycled.

Recycling the packaging reduces the use of raw materials in the manufacturing process and also reduces the amount of waste in land-fill sites.

Optional components

- A 314, mesh tray for reprocessing lumened instruments, with lockable lid, holding lugs and 1 injector manifold with Luer Lock ♀ connector for Luer Lock ♂, height 47 mm, width 151 mm, depth 225 mm
- E 197, mesh tray with lockable lid and 2 holding lugs for instruments, height 42 mm, width 150 mm, depth 225 mm
- E 447, Luer Lock Adapter ♀ for Luer Lock ♂

Further components are available from Miele as optional extras.

en - Warnings and safety notes

To avoid the risk of accidents and damage to this mobile unit please read these instructions carefully before using it for the first time. Keep these instructions in a safe place where they are accessible to users at all times.

⚠ Please also read the operating instructions for your dishwasher or the machine and pay particular attention to the Warning and Safety instructions.

- ▶ The mobile unit is approved solely for the applications specified in the “Areas of application” chapter of these operating instructions. Applications that go beyond those listed and use of the mobile unit with additional components from the Miele product range are described in the operating instructions for the components or are to be agreed with Miele.
- ▶ New load carriers must be cleaned in the washer-disinfector without a load prior to first use.
- ▶ Inspect all mobile units, baskets, modules and inserts daily as described in the “Maintenance” section in the Operating instructions for the washer-disinfector.
- ▶ For your own safety, wear gloves and, if necessary, a surgical mask and protective goggles when sorting the instruments and unloading the machine. To avoid recontamination, wear clean gloves when unloading the machine.
- ▶ Ophthalmic instruments should only be reprocessed in a machine **specifically dedicated for this application**. This way there is no risk of residues from other disciplines being deposited in narrow lumened hollow instruments.
Synthetic cover nets, e.g. Miele cover nets A 2 or A 3, must not be used in this machine.
- ▶ When reprocessing ophthalmic instruments do not dispense rinsing agent.
- ▶ The reprocessing programme used should include thermal disinfection, and the process chemicals used must be suitable for the material quality of the load items.
- ▶ The last but one and the final rinse should be carried out with fully demineralised water wherever possible.

⚠ Ophthalmic instruments should only be reprocessed in a washer-disinfector **specifically dedicated for this application**. This way there is no risk of residues from other disciplines being deposited in narrow lumened hollow instruments.
Synthetic cover nets, e.g. Miele A2 or A3 cover nets, must not be used in this washer-disinfector.

⚠ When reprocessing ophthalmic instruments do not dispense rinsing agent.

Programme selection

The A 207 mobile injector unit may **only be used in the** OphthaTrays A207 **programme**. This programme uses thermal disinfection.

To check before loading the machine and before starting a programme

- Are the silicone hoses securely screwed on?

⚠ Make sure that all screw connectors are fitted with silicone hoses with Luer-Lock connectors or blind stoppers to ensure that all fittings in use are supplied with sufficient standardised pressure. Damaged silicone hoses and Luer-Lock connectors must not be used.

Unused silicone hoses with Luer-Lock connectors must be secured in the perforated plate or be replaced with blind stoppers.

- Are the spray arms able to rotate freely?
- Check that the mobile unit is correctly docked to the water path in the chamber.

Loading ophthalmic instruments

For validation purposes please make sure that you follow the loading instructions given on the template.

⚠ Risk of injury caused by load.

There is a risk of injury when loading and unloading due to possible sharp edges, rims or pointed ends.

To minimise the risk of injury, loading should take place from the back forwards and unloading should take place in the reverse order.

⚠ The wash water used for reprocessing narrow lumen instruments needs to be filtered. The injector manifold must be fitted with tubular filter A 800. The tubular filter is supplied with its own fitting instructions.

⚠ Before connecting hollow instruments, please make sure that all hollow sections and cannulae are free of blockages. Otherwise the cleaning solution will not be able to access all areas properly. The cannulae and other hollow instruments should be flushed through with a spray of water as soon as possible after use on patients, preferably whilst still in the theatre. This is to prevent residual blood or medication causing a blockage in the narrow lumen and endangering internal cleaning.

⚠ Risk of infection due to insufficient disinfection.

Load items that become detached from the irrigation connector during reprocessing will not be sufficiently disinfected on the inside. Any items that become detached from the irrigation connector during reprocessing must be reprocessed again.

⚠ Instruments must not touch each other during reprocessing.


Suitable mesh trays

Non-Miele mesh trays, cassettes and other containers can be used for reprocessing instruments in Miele washer-disinfectors if they have been risk assessed and specified as suitable for such use and if reprocessing procedures have been successfully validated with their use.

Loading lumened instruments

- Arrange the instruments in the mesh trays according to their manufacturer's instructions.


When loading the instruments, please note:

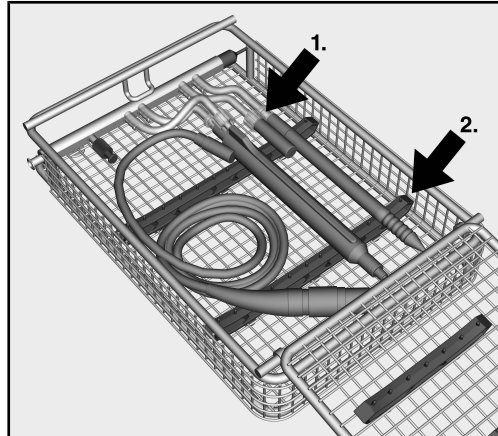The example shows the A 314 mesh tray

1. When processing hollow instruments with multiple connectors (rinsing and suction hand pieces, for example), each connector needs to be attached to a hose connector in the mesh tray. Hoses must not be twisted or kinked. If hoses are twisted or kinked, less wash water will flow through the instruments and this will affect reprocessing results.
2. Lightweight and delicate instruments must be secured in silicone or plastic holders. Space the holders in the mesh tray according to the size of the instruments.

Loading mesh trays without injector manifolds

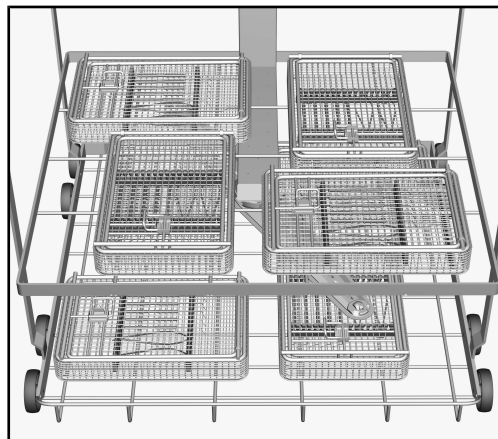- Arrange the mesh trays centrally in the lower levels of the A 207 mobile unit as shown.

The mesh trays must not be stacked on top of each other.

en - Areas of application


Loading mesh trays with injector manifolds

- Load mesh trays with injector manifolds upright in the mesh tray carrier A 311. The water connections must point towards the injector manifold on the mobile unit. Load from the back to the front, and unload from the front to the back.

The mesh trays must not be placed loosely in the mobile unit or stacked on top of each other.

- Connect the mesh trays to the silicone hoses on the A 207 mobile unit (the A 314 mesh tray is shown in the example here).

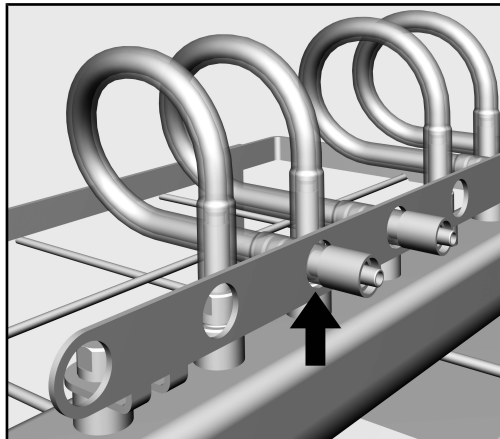- Secure the ends of any unused hoses in the perforated plate.

⚠ Otherwise unused silicone hoses will be spun around by the force of the water. This could cause damage to the load or to the machine.

Unused hoses must not be sealed with stoppers.

A 314 mesh tray

Always check the following before and after using the A 314 mesh tray:

- Is the injector manifold screwed tight to the mesh tray?
- Are the hose connections at the injector manifolds free of blockages?
- Are the silicone hoses on the mesh tray connected securely to the injector manifold?
- Are there any signs of wear on the silicone hoses or the Luer Lock adapters?

Replace the silicone hoses on the mesh tray after approx. 25 programme sequences or as soon as there are visible signs of wear.

Testing lumened instruments within the framework of quality assurance tests

 Risk of chemical burns from alkaline residues.
Even minimal amounts of alkaline detergent residues can cause chemical burns to the retina.
After machine reprocessing, check the pH value of the residual water in the instruments.

- Residual water must be expelled from lumened instruments using medical compressed air. This procedure also verifies that lumens are clear.
- Check the pH value of the water droplets that come out of the instruments using indicator paper with increments of ≤ 0.5 . The pH value should be between 5 and 8.

If the test result is negative


- In the event of a blockage or reduced clearance, remove the offending particles or blockage manually.

Then, even if the pH value deviates from neutral:

- Reprocess the lumened instruments.

If the pH is still not neutral after reprocessing, contact Miele to arrange for further analysis of the problem.

Cleaning the A 800 tubular filter

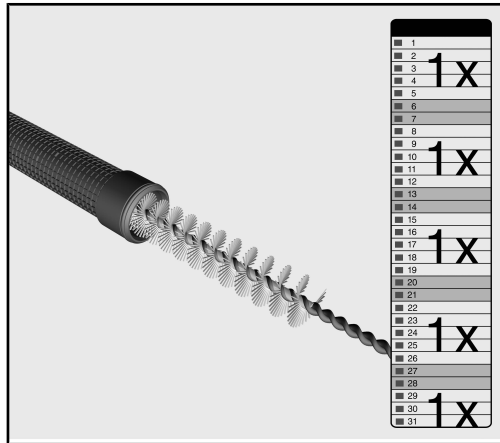
The tubular filter must be cleaned **at least once a week**.

We recommend cleaning the tubular filter after every 10 programme sequences.

- Follow the cleaning instructions in the operating instructions for the tubular filter.

Test point for measuring water pressure

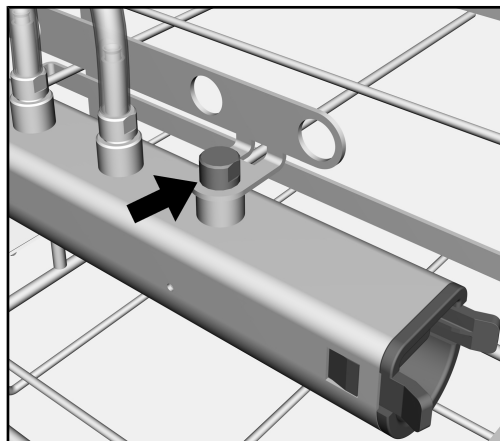The connection at the front is used as the test point for measuring water pressure. This connection point can be used for testing pressure within the framework of performance and validation testing in accordance with EN ISO 15883.

- To measure water pressure replace the blind stopper with a Luer-Lock adapter, e.g. E 447.

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